


Innova-YSQ Yellow Square Edge


Now Filters both UV & Blue Light

- 360° Square Edge Design
- Yellow Cromophore
- Singlepiece IOLs

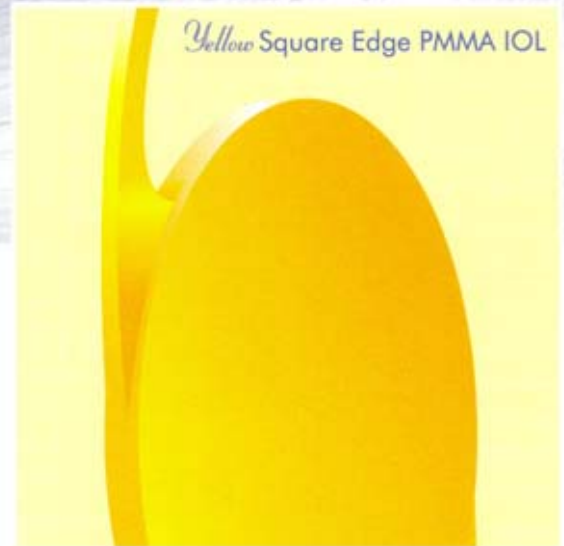


The Innova Yellow PMMA Lenses are having ability to filtering blue light and to provide protection that could reduce the risk of AMD among cataract patients. As you know that our natural lens has ability to filter blue light comes from UV but artificial lens doesn't. The Innova Yellow Cromophore PMMA Lens filter the blue light & protect the retina from the same.

Technical Specification

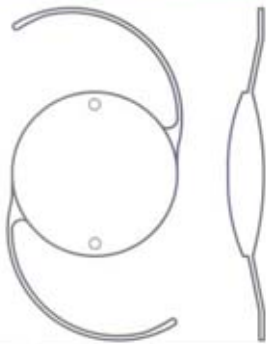

Material	YELLOW PMMA
Manufacture	Lathe Cut, Tumble Polished
Model No.	YSQ 600
Optic Style	Equi Convex
Optic Size	6.00 mm Diameter
Diametral	12.50 mm Diameter
Haptic Style	Modified C-Loop,
Position Hole	2
A "Constant"	118.0
Diopters	8.00 to 30.00

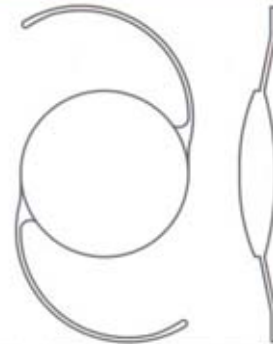

Material	YELLOW PMMA
Manufacture	Lathe Cut, Tumble Polished
Model No.	YSQ 525
Optic Style	Equi Convex
Optic Size	5.25 mm Diameter
Diametral	12.00 mm Diameter
Haptic Style	Modified C-Loop,
Position Hole	0
A "Constant"	118.2
Diopters	8.00 to 30.00



OMNI LENS PVT. LTD.

Office : 5, "Samrudhi", Opp. Sakar-III,
Navrangpura, Ahmedabad-380 014. (India)
Tele.: 27544013, 30921703 Fax : 0091-79-27544035
E-mail : bgpgroup1@icenet.co.in
Web: www.bgpgroup.com